

BestTM Cyclotron Systems



Best Cyclotron Systems provides 15/20/25/30/35/70 MeV Proton Cyclotrons as well as 35 and 70 MeV Multi-Particle (Alpha, Deuteron & Proton) Cyclotrons

- Currents from 100uA to 1000uA (or higher) depending on the particle beam are available on all BCS cyclotrons
- Best 20u to 25 and 30u to 35 are fully upgradeable on site

| Cyclotron | Energy (MeV) | Isotopes Produced |
|------------------------|--------------|--|
| Best 15 | 15 | ¹⁸ F, ^{99m} Tc, ¹¹ C, ¹³ N, ¹⁵ O, ⁶⁴ Cu, ⁶⁷ Ga, ¹²⁴ I, ¹⁰³ Pd |
| Best 20u/25 | 20, 25-15 | Best 15 + ¹²³ I, ¹¹¹ In, ⁶⁸ Ge/ ⁶⁸ Ga |
| Best 30u (Upgradeable) | 30 | Best 15 + ¹²³ I, ¹¹¹ In, ⁶⁸ Ge/ ⁶⁸ Ga |
| Best 35 | 35-15 | Greater production of Best 15, 20u/25 isotopes plus ²⁰¹ Tl, ⁸¹ Rb/ ⁸¹ Kr |
| Best 70 | 70-35 | ⁸² Sr/ ⁸² Rb, ¹²³ I, ⁶⁷ Cu, ⁸¹ Kr + research |



Assembly of a Best 35 MeV Cyclotron at Best Theratronics facility, Ottawa, Ontario, CA


Installation of Best 70 MeV Cyclotron at Italian National Laboratories (INFN), Legnaro, IT



BestTM Particle Therapy

400 MeV Rapid Cycling Medical Synchrotron for Proton-to-Carbon Heavy Ion Therapy:

- Intrinsically small beams facilitating beam delivery with precision
- Small beam sizes – small magnets, light gantries – smaller footprint
- Highly efficient single turn extraction
- Efficient extraction – less shielding
- Flexibility – heavy ion beam therapy (protons and/or carbon), beam delivery modalities

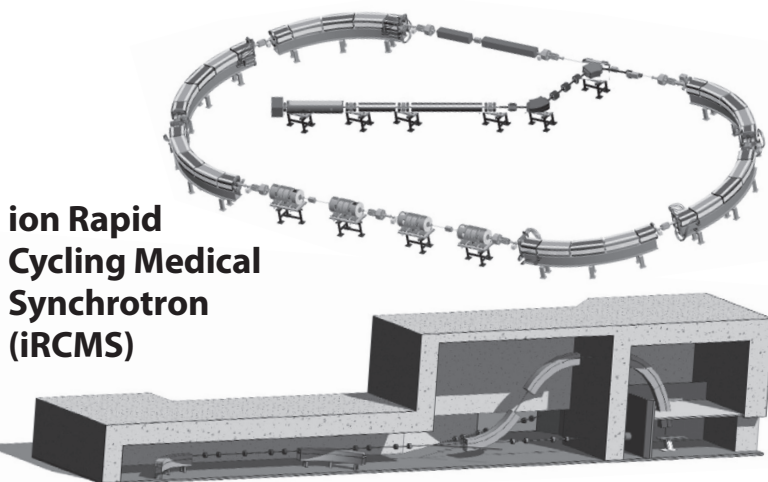

ion Rapid Cycling Medical Synchrotron (iRCMS)

TeamBest Companies © 2017–2018

* Best iRCMS is under development and not available for sale currently.

www.bestcyclotron.com • www.bestproton.com • www.teambest.com

USA tel: 703 451 2378 800 336 4970 • CANADA tel: 604 681 3327 866 909 4647